

signia

Life sounds brilliant.

REPLICATING NATURE.

Signia Nx with OVP™ for the most natural own voice
and highest acceptance



Products

No Software update

- ClickSleeves 2.0

Availability

March 2018

Products

Software CXX8.5

- Pure Charge & Go Nx
- Motion 13 P Nx
- Pure 10 Nx
- Pure 13 Nx Update [T]
- Pure CROS Nx
- Streamline Mic

Availability

April 2018

Products

Software CXX8.5.5

- Insio CIC/IIC Nx
- Silk Nx (not OHS approved)
- Silk CROS Nx

Availability

June 2018

Signia Nx Products

StreamLine Mic

StreamLine Mic enables direct audio streaming from any Bluetooth device

The new StreamLine Mic can be connected to virtually any Bluetooth-enabled smartphone and establishes various direct audio streaming options:

- Audio streaming from Android smartphones
- Provides full hands-free phone calls
- Enjoy music in stereo sound
- Hear in challenging situations with remote microphone mode



Compatibility:

- Pure 312 7/5/3 Nx
- Pure 13 7/5/3 Nx
- Pure C&G 7/5/3 Nx
- Motion 13 7/5/3 Nx

Signia Nx Products

StreamLine Mic - Overview

StreamLine Mic offers intuitive onboard control and status elements


Multi function button

- Switch on (press 3 sec) and off (press 6 sec)
- Toggle between two modes:
 - Bluetooth Phone Mode
 - Remote Microphone Mode

Static directional microphones

VC / Program up VC / Program down

- integrated Microphone for
 - Handsfree phone calls
 - Remote microphone functionality
- Accept/decline calls
- Volume increase
- Volume decrease

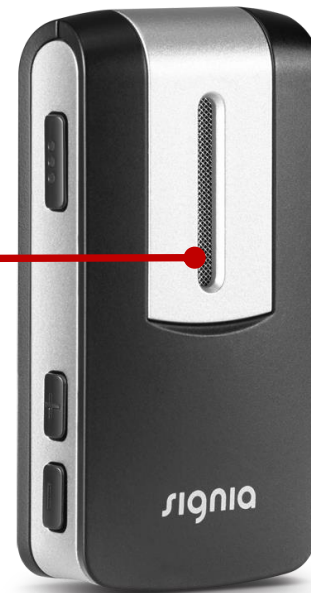

Power indicator

Status indicator



Signia Nx Products

Pure Charge&Go Nx (April 2018)

Pure Charge&Go Nx

The ultimate RIC solution which combines all the benefits of the recently launched Pure 312 Nx (OVP, Ultra HD e2e, Bluetooth, Telecare 3.0) with our successful Li-Ion technology. **Available in the 7nx/5nx/3nx performance levels.**

Same Li-Ion module and inductive charger as provided by Cellion primax
The overall size comparable in size to Pure 312 Nx.



Feature-benefit-value matrix

Feature	Benefit	Value
Inductive rechargeability with integrated design	No purchasing and handling of disposable batteries. Robust devices with less moving parts and openings.	Just charge and go with convenient handling and reliable hearing.
Bluetooth connectivity with Li-ion power	High capacity cell provides extensive Bluetooth direct streaming.	Convenient and carefree enjoyment of TV, music and phone calls.

Signia Nx Products

Motion 13 P Nx (April 2018)

Motion 13 P Nx

Motion 13P Nx extends the fitting range of the Signia Nx platform to severe to profound hearing losses with the convenience of a BTE. **Available in the 7nx/5nx/3nx performance levels.**

- Similar design (look and feel) compared to Motion 13 Nx (same T-coil battery doors and coupling system)
- Improved fitting range / matrix of power BTE - more high frequency amplification compared to Motion P (Matrix: 135/77dB (prel) compared to 131/70dB for Motion P px)



Feature-benefit-value matrix

Feature	Benefit	Value
Increased high frequency amplification	Maximizes essential speech understanding for greater hearing losses.	Enjoy conversations with clearer speech
Bluetooth and T-coil	Improves speech understanding through connectivity	Enjoy clearer meetings, TV, music and phone calls.

Signia Nx Products

Pure 10 Nx (April 2018)

Pure 10 Nx

Pure 10 Nx is targeted at clients with very high cosmetic needs. Pure 10 Nx compromises on functionality (i.e. no connectivity) to achieve the highest discreetness possible. **Available in the 7nx/5nx/3nx performance levels.**

Pure 10 is the successor of Ace primax and remains w/o HD e2e or Bluetooth.

- Features improved sound quality due to Signia Nx sound processing.
- Smooth flat housing w/o the pb function (optional PU available)



Feature-benefit-value matrix

Feature	Benefit	Value
Non-wireless with no on-board controls	Reduced connectivity limits user errors and interaction.	A simple device that you don't have to adjust.
Ultra-discreet design with TwinMic	Smallest behind-the-ear device with directionality	A hidden device with high performance in noise.

Signia Nx Products

Pure 13 Nx Update (April 2018)

Pure 13 Nx Update

The Pure 13 Nx will be made compatible to the T-Coil Battery Door. **Available in the 7nx/5nx/3nx performance levels.**

This is a new product with updated hardware (contacts), thus new part/order numbers are required! The existing Pure 13 Nx cannot be upgraded and will be phased out with the launch of the updated product.



Feature-benefit-value matrix

Feature	Benefit	Value
Longest wearing time in its class	Most energy-efficient streaming thanks to the high-capacity 13-size battery	The wearer can enjoy connectivity and streaming without any concern on battery life

Signia Nx Products

Cros Pure Nx (April 2018)

Cros Pure Nx (312)

Wireless CROS/ bi-CROS transmitter for use with the Signia Nx portfolio. Binaural processing for CROS Pure (not available before in primax) to help understanding in noise for patients who really need it (big problem when having single-sided-deafness and certain hearing loss on the good side)

- Bluetooth connectivity for CROS Pure. (World's first MFI CROS solution)
- Same housing as Pure 312 Nx and compatible with all Bluetooth Nx products (Pure 312 Nx, Pure 13 Nx, Motion 13 Nx, Motion 13P Nx)



Feature-benefit-value matrix

Feature	Benefit	Value
Bluetooth and CROS	CROS fittings are compatible with Bluetooth enabled hearing aids.	Enjoy music and phone calls direct from smart phones.
Ultra HD e2e	Narrow directionality improves speech understanding and eWindscreen binaural reduces unwanted wind noise.	Easier conversations even in noisy situations and more listening comfort in wind.

Silk Nx (not OHS endorsed)

For wearers wanting to keep their hearing loss hidden to everyone around them, Silk Nx offers the most discreet, instant-fit solution for the most natural hearing. **Available in the 7nx/5nx/3nx performance levels.**

- Instant-fit, discreet (size 10) CIC solution (20% smaller than primax)
- Click sleeve 2.0 compatibility (incl new XS size providing higher fit rate)
- Ultra HD e2e, touchControl App, TeleCare 2.0
- Does not feature Bluetooth, does not feature OVP



Feature-benefit-value matrix

Feature	Benefit	Value
20% smaller size with darker face-plates	Increased fitting rate and further decreased visibility.	A discreet hearing solution no one will see
Instant fit design with click Sleeves	Ready-to-wear and secure fit for almost every ear.	The most discreet hearing solution immediately and hassle-free.

Signia Nx Products

Insio CIC / IIC Nx (June 2018)

Insio CIC / IIC Nx

Nx custom IIC / CIC introduces custom-made solutions completely in the canal with the benefits of the Nx sound engine w/o Bluetooth connectivity. User benefit from binaural connectivity features such as binaural One Mic for better speech understanding in noise. **Available in the 7nx/5nx/3nx performance levels.**

- IIC size: New small size faceplate (“I”) & thin shell material process & smaller microphone
- IIC & CIC (T5-FP) size: smaller components, wireless coil & chip
- Full Nx binaural features



Feature-benefit-value matrix

Feature	Benefit	Value
Smaller components	Smaller instruments and an increased fit rate	A discreet and comfortable hearing aid no one will see.
Binaural OneMic directionality	Provides better speech understanding in noise.	Enjoy conversations in noise with clearer speech

Signia Nx Products

Cros Silk Nx (June 2018)

Cros Silk Nx (10)

Wireless CROS/ bi-CROS transmitter for use with the Signia Nx portfolio. This is the world's first CIC solution for single sided deafness offering the most discreet natural experience for CROS wearers




- CROS Silk does not feature OVP or Bluetooth functionality
- CROS Silk is compatible with the following products: Silk Nx, Insio CIC Nx and Insio IIC Nx



Feature-benefit-value matrix

Feature	Benefit	Value
CIC CROS	CROS functionality for completely in the canal devices	A discreet hearing solution no one will see

Signia Nx Products Technology Summary

	RIC 10	RIC 312	RIC 13	RIC LI	BTE 13	BTE 13 P	Custom	Click ITE	CROS	
7/5/3	Pure 10	Pure 312	Pure 13	Pure Charge& Go	Motion 13	Motion 13P	CIC/IIC	Silk CIC	CROS Pure 312	CROS Silk
	-	✓	✓	✓	✓	✓	✓	✓	✓	✓
	-	✓	✓	✓	✓	✓	-	-	-	-
	-	✓	✓	✓	✓	✓	-	-	✓	-

Silk will not be OHS approved

REPLICATING NATURE.

Signia Nx with OVP™ for the most natural own voice
and highest acceptance

